



AESTHETICS

Non-surgical innovations

ADVANCED TECHNOLOGY, AFFORDABLE PRICES, NO COMPROMISE ON SPECIFICATION

INTRODUCING




Dermaforce

ADVANCED RADIOFREQUENCY MICRONEEDLING



THE 3D APPROACH

A multi layered approach to facial rejuvenation

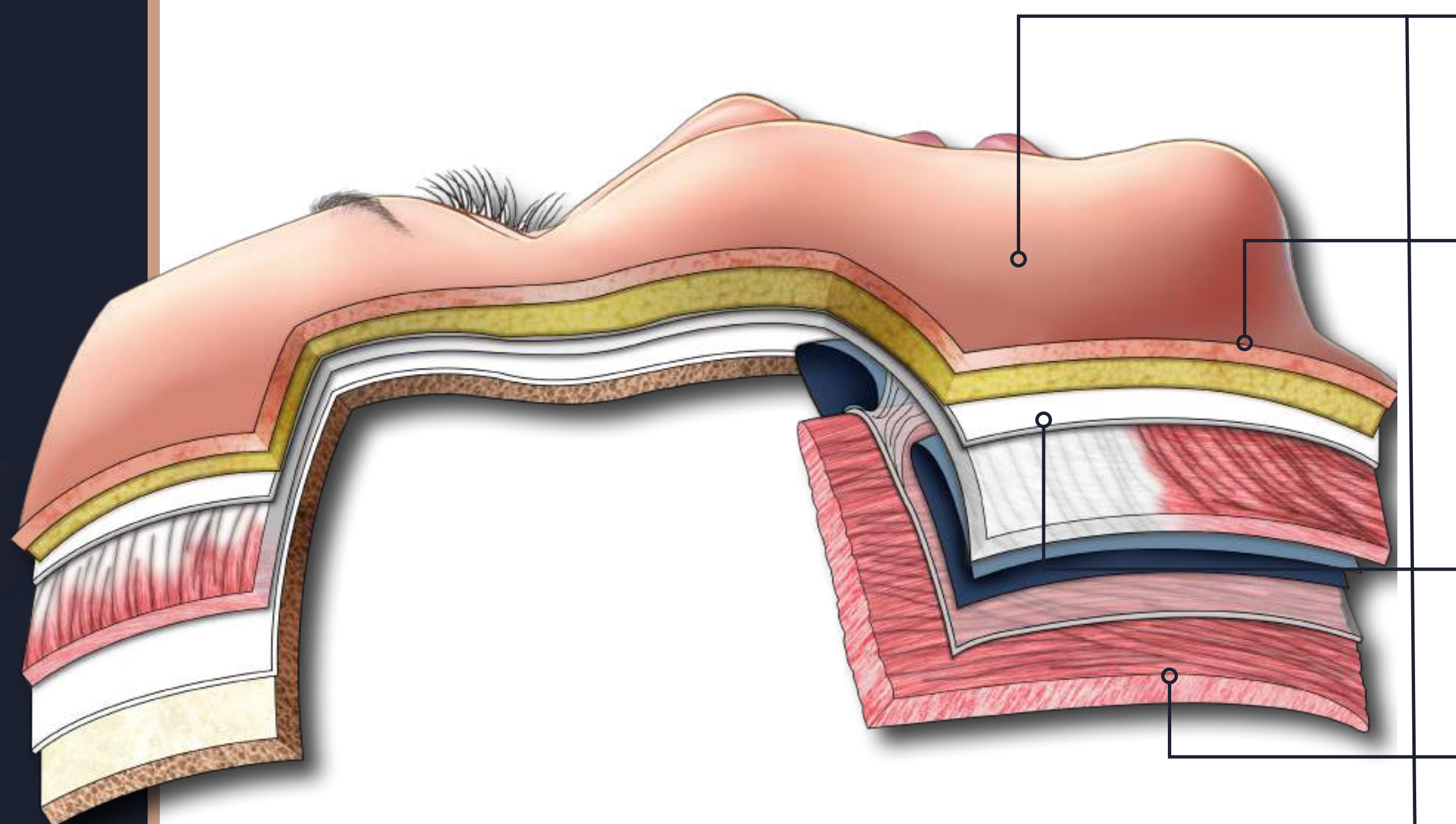

1. EPIDERMIS

- **HYDRO₂ FACIAL NON R/F:** Overall skin health
- **ABLATIVE LASERS:** pigmentation/skin resurfacing
- **VJUVE MD:** CO₂ Laser skin resurfacing

2. DERMIS

- **3D DERMAFORCE® / TRIMAX MD:** Microneedling RF Collagen Stimulation
- **3D HYDRO₂ / 3D-LIPO ULTIMATEPRO:** Radio Frequency Skin tightening
- **VJUVE MD- CO₂ LASER:** Skin resurfacing

3. SMAS LAYER UPWARDS

- **3D HIFU:** High Intensity focused Ultrasound non-surgical facelift

4. SMAS LAYER UPWARDS

- **EMLIFT MD:** technology accelerated muscle stimulation

5. SMAS LAYER UPWARDS

- **LED:** Regeneration

WHY CHOOSE 3D DERMAFORCE® FOR YOUR BUSINESS?

What makes the 3D Dermaforce® Different?

- Specifications more advanced than market-leading devices
- at a fraction of the cost!
- Advanced, Unique R-Flow technology - The safest device on the market, controlled by vacuum delivery engineered to eliminate risk of cross contamination
- Achieve Fantastic results in just 3 x 30 Minute sessions!
- High Hourly Revenue
- Multiple Treatment Options from just one device!
- Cross treatment protocol opportunities
- As featured in the press and media



INTRODUCING THE 3D DERMAFORCE®




THE TECHNOLOGY

A highly effective and versatile method of skin rejuvenation

3D Dermaforce® Microneedling/RF is an advanced method of skin rejuvenation that stimulates the body's natural healing process to encourage collagen production in the dermis resulting in rejuvenated and tightened skin.

The Power of Advanced RF Microneedling RF is making waves in the industry - proven to reduce fine lines, wrinkles, and lift the skin, this minimally invasive treatment also stimulates the production of collagen and elastin for firmer, more youthful skin.

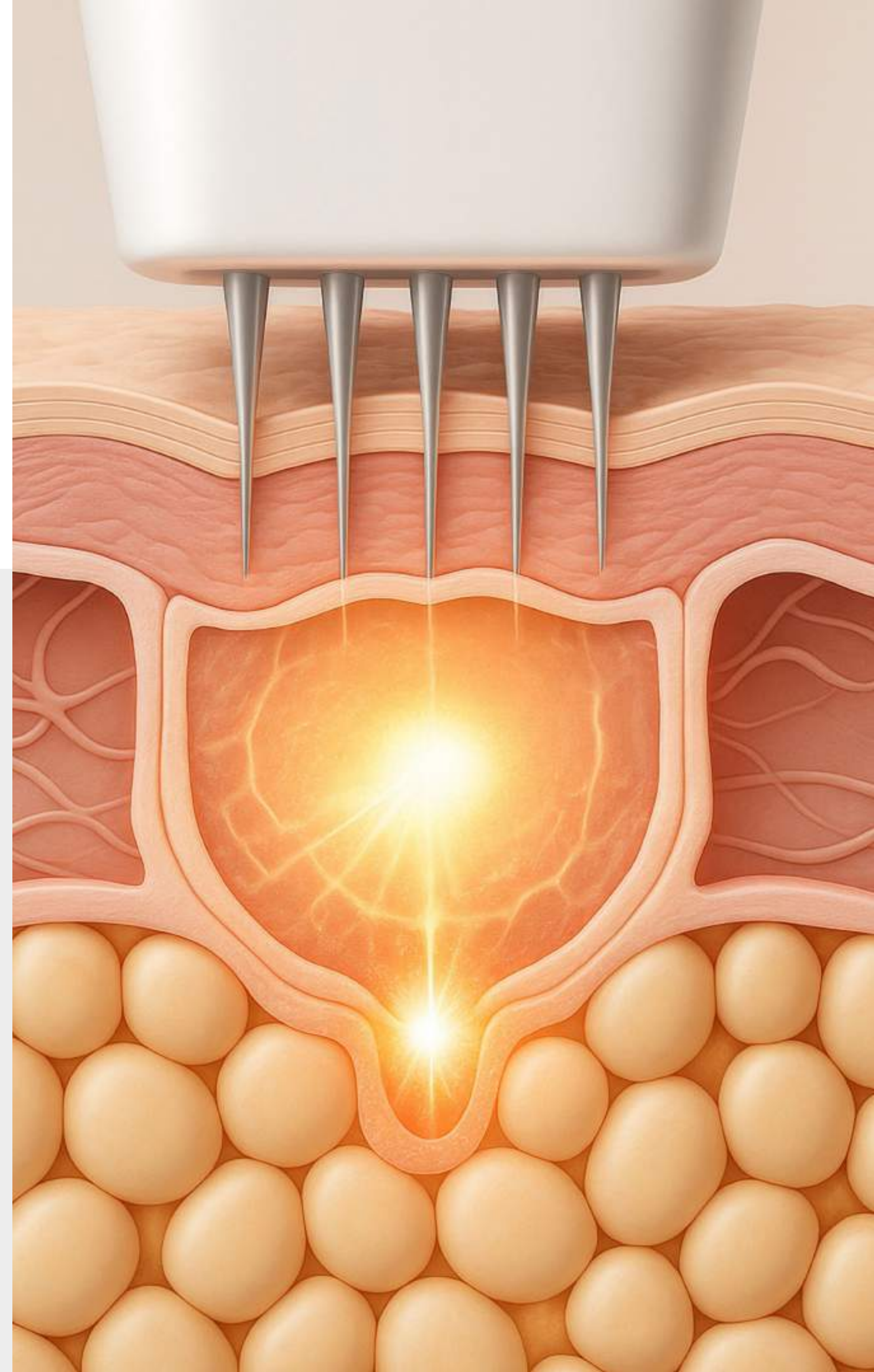


THE TECHNOLOGY

A Highly Effective and Versatile Method of Skin Rejuvenation

3D Dermaforce® Microneedling/RF is an advanced method of skin rejuvenation that stimulates the body's natural healing process to encourage collagen production in the dermis resulting in rejuvenated and tightened skin.

The Power of Advanced RF Microneedling RF is making waves in the industry - clinically proven to reduce fine lines, wrinkles, and lift the skin, this minimally invasive treatment also stimulates the production of collagen and elastin for firmer, more youthful skin.

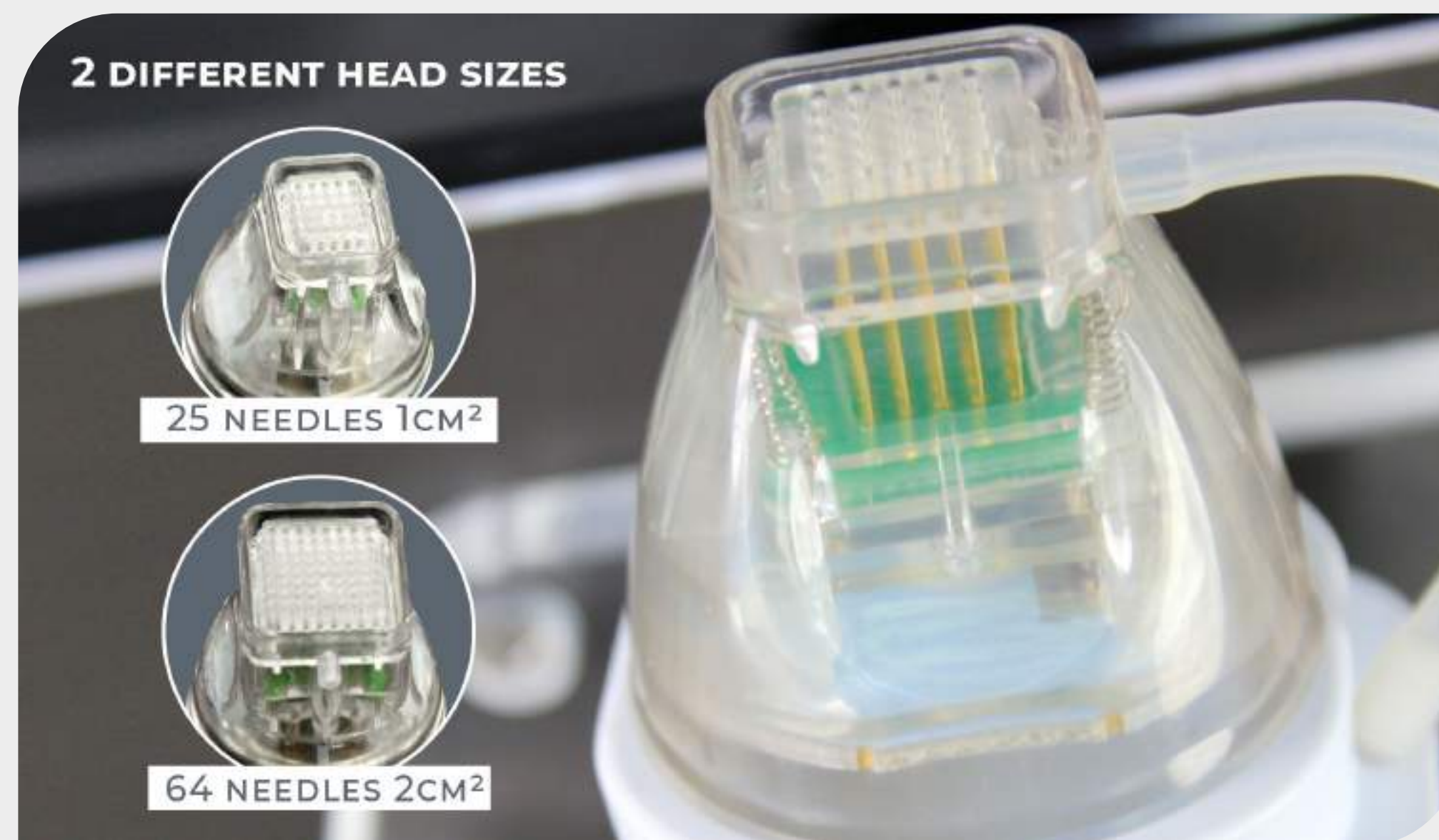


THE TECHNOLOGY

A Highly Effective and Versatile Method of Skin Rejuvenation

During treatment, the radio frequency energy is delivered through the insulated needles allowing precise delivery to the dermis while protecting the epidermis layers of the skin.

The radio frequency wave penetrates the basal melanocyte barrier, heating the collagen fibres in the dermis to 55°C-65°C.

This stimulates the body's natural healing process, increasing cell turnover as the body heals – resulting in the production of collagen and elastin in these areas.

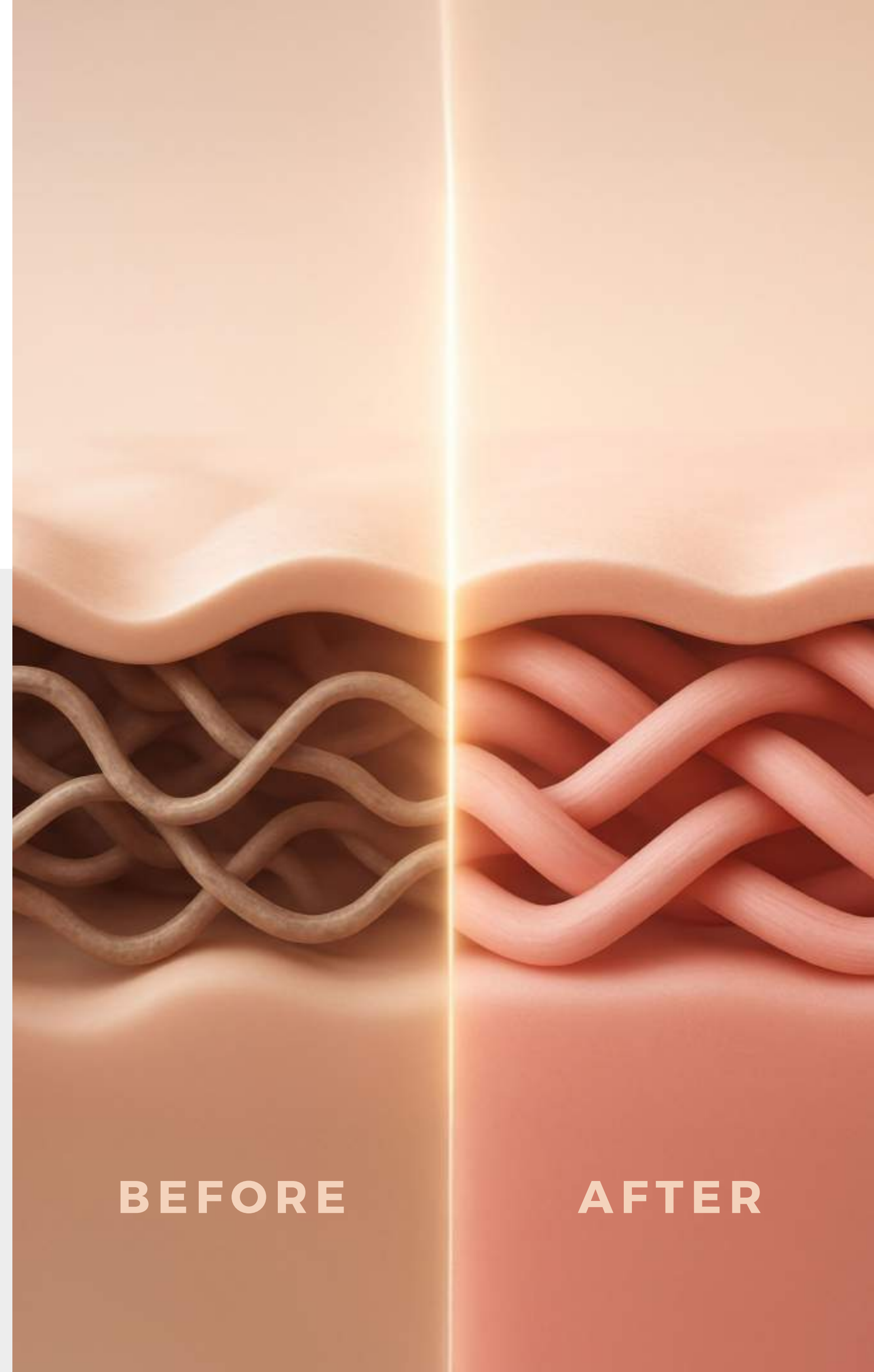

BEFORE

AFTER

INTRODUCING R-FLOW® TECHNOLOGY



Introducing 3D Dermaforce® with Patented R-Flow Technology eliminating any risk of cross contamination.

The 3D Dermaforce® device is setting new standards in reliability with its exclusive, patented R-Flow® Technology – a breakthrough developed by our in-house engineering team.

The only device in the market with engineered R-Flow Technology – developed and installed by our in-house engineering team, Patent pending, exclusive to 3D. Council by council, restrictions are being enforced to ensure that RF Microneedling devices comply with the highest level of safety concerns – so far, there are only a handful of council which are enforcing this, but the number is growing.

This advanced engineering development means that 3D Dermaforce® is the only RF Microneedling device on the market powered by controlled vacuum delivery to eliminate any risk of cross contamination – achieved by superior seal systems, improved filtration and reverse flow technology.

INTRODUCING R-FLOW® TECHNOLOGY

*Introducing 3D Dermaforce® with Patented R-Flow Technology
eliminating any risk of cross contamination.*



TREATMENT VERSATILITY

This Treatment is so versatile that it can be used to target a number of concerns on the Face and Body, including:

- ✓ Stretch marks
- ✓ Lines and wrinkles
- ✓ Overall skin texture
- ✓ Pigmentation
- ✓ Sun damaged skin
- ✓ Pitting of skin
- ✓ Sagging skin
- ✓ Enlarged pores
- ✓ Scarring



Loose and
sagging skin




Fine lines
and wrinkles



Sun damage/
Pigmentation



Stretch
Marks

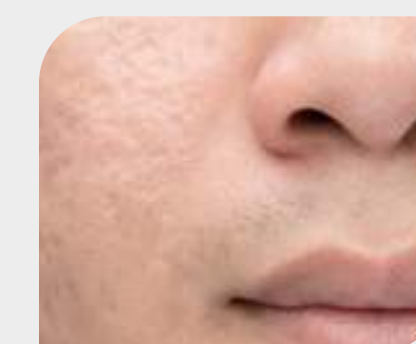

Pitted skin/
skin texture

TREATMENT VERSATILITY

Popular Treatment Areas Include:

FOR THE FACE:

- Crows Feet / Eyes
- Nasolabial fold
- Jowls and Jawline
- Marionette lines
- Neck and Decollete
- Necklace Lines
- Forehead

FOR THE BODY:

- Abdomen
 - Mummy Tummy
 - Surgical scar revision
 - Skin Tightening
- Knees
- Hands
- Bingo wings / upper arms
- Inner thighs
- Elbows



GREAT RESULTS



Photo credit: Skin You Will Love

GREAT RESULTS



GREAT RESULTS



GREAT RESULTS




GREAT RESULTS



Photo credit: This is Me Aesthetics

GREAT RESULTS

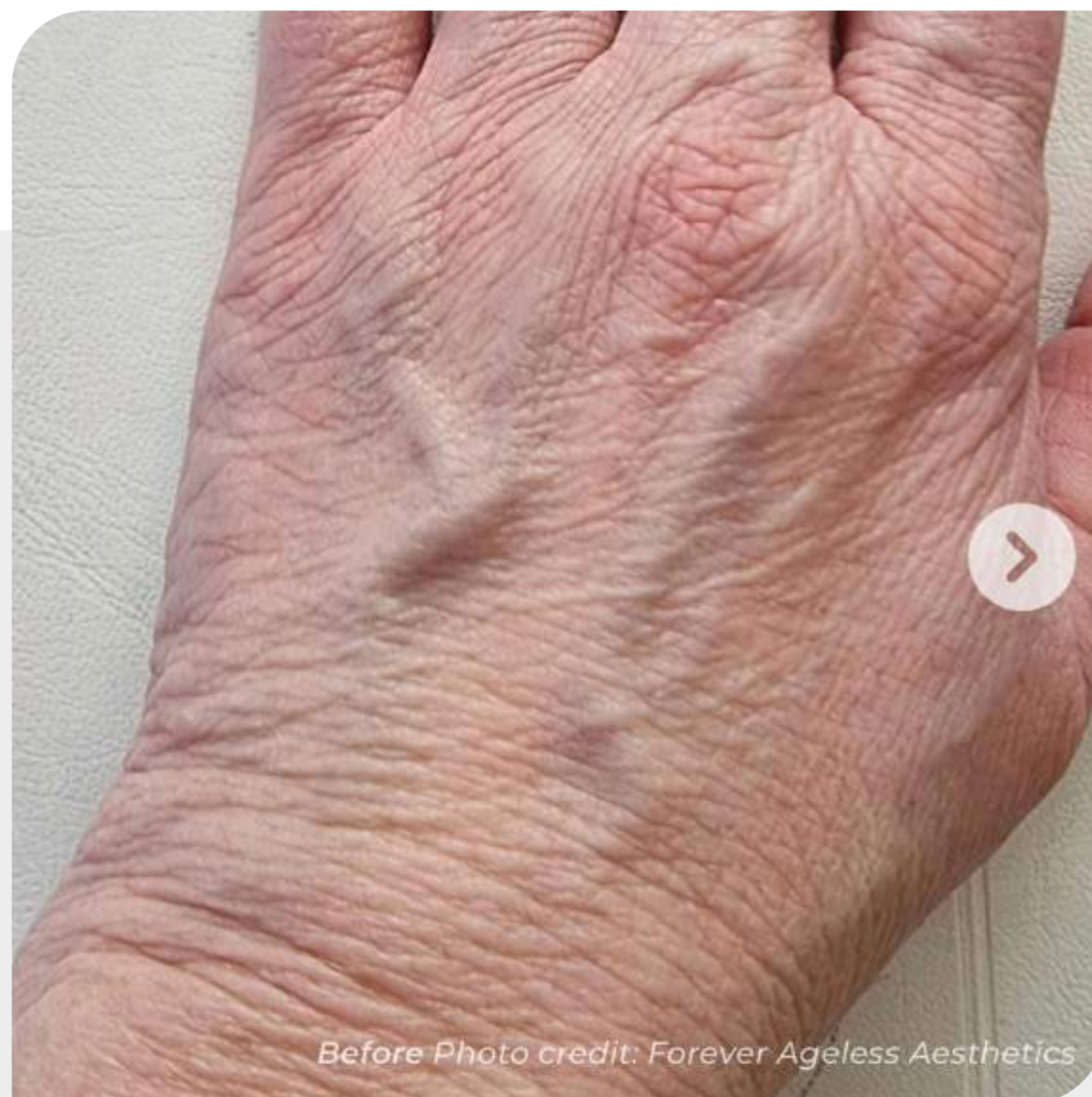

3D DERMAFORCE®

KEY FEATURES:

- Insulated needles for precise delivery to target the dermis and protect the epidermis.
- Accurate, controlled delivery.
- Stepping motor for patient comfort, allowing needles to be smoothly inserted into the skin.
- Precise Depth Control 0.5mm–3.5mm (0.5mm step).
- Fine needles of 0.12mm for improved client comfort.
- Vacuum feature for auto-injection offers a more comfortable treatment.
- Treatment head size – 36 needles, 64 needles.
- Gold plated needles, suited to clients with metal allergies.
- Ergonomic handpiece design, including trigger button.

THE TREATMENTS



FACE – BROW LIFTING



FACE – FINE LINES AND WRINKLES




BODY – STRETCH MARKS




BODY – SKIN TIGHTENING

LOVED BY PRESS AND CELEBRITIES

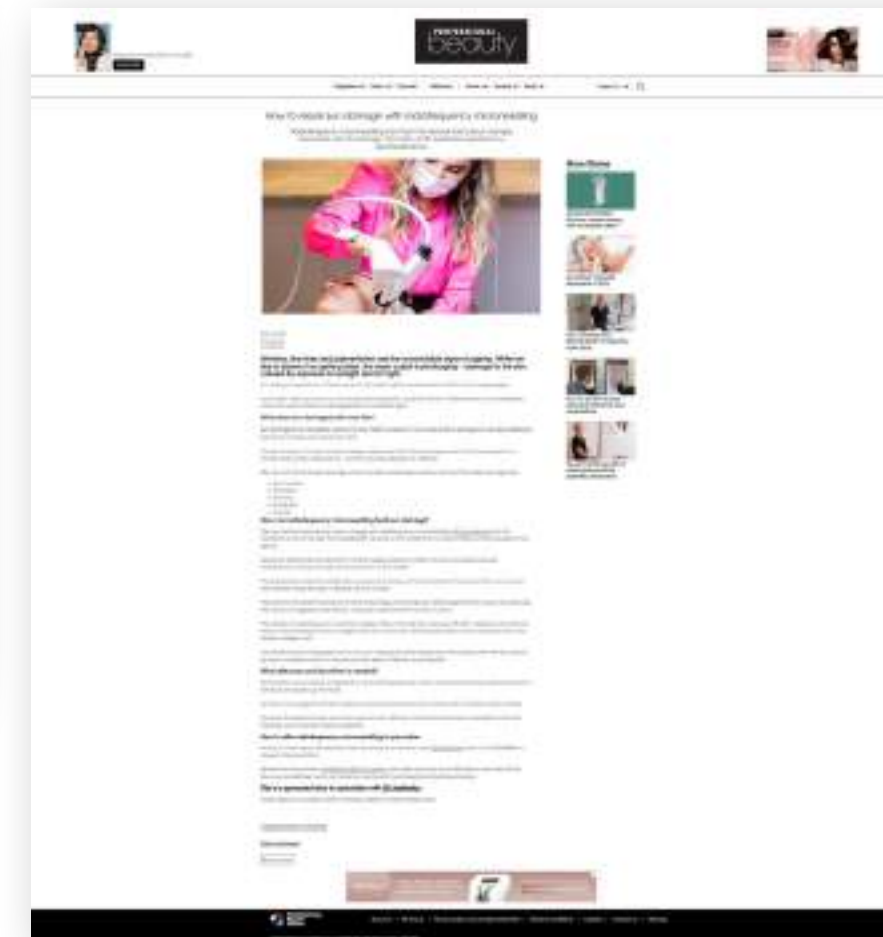

GREAT TESTIMONIALS

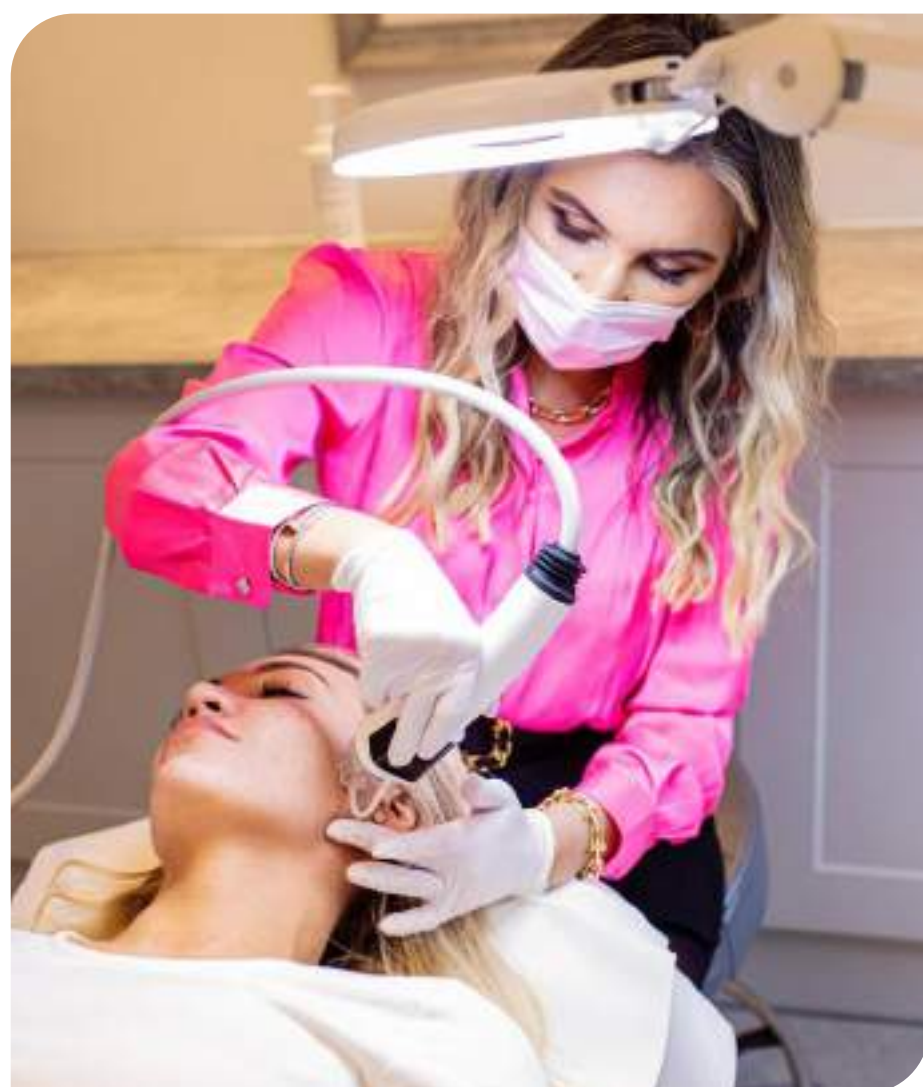

“

As their Key Option Leader, I am really proud to partner with 3D Aesthetics. I think that they are innovative, are constantly moving forward, and I have really enjoyed working with their products and devices within my own practice.

”

Dr Jonquille Chantrey



“

We partnered with 3D Aesthetics and Dermaforce® to expand our treatment offerings, and it's been a game-changer. Not only did we generate £30,000 in revenue in just 2 weeks, but the support and training they provided made it all so easy. This partnership has transformed our clinic!

”

Alison Telfer, Glass House Clinic



“

We at the nurses network choose to partner with 3D Aesthetics because they are trusted, credible suppliers. “Their post-purchase support is incredible, and their results are un-matched.

”

Megan Hope, MCH Clinic

TREATMENT PRICING

TREATMENT TYPE	Cost of Treatment	Recommended Price Per Treatment	Recommended Price Per Course of 8
Full Face Treatment (30 mins)	£30	£250-400	£750-1200
Full Face & Neck Treatment (45 mins)	£30	£300-500	£900-£1500
Crows Feet / Eye Treatment (15 mins)	£30	£150-200	£450-600
Tummy Tightening Treatment (45 mins)	£30	£300-500	£900-£1500
Bingo Wings / Other Skin Tightening (30 min)	£30	£250-400	£750-1200

RECOMMENDED RETURN ON INVESTMENT

Treatment Revenue Potential

TREATMENT	Treatment Price	Weekly Treatments	Monthly Revenue	Month 3	Month 6	Month 9
Full Face Treatment (Single Treatment)	£400	1	£1,600	£4,800	£9,600	£14,400
	£400	2	£3,200	£9,600	£19,200	£28,800
	£400	5	£8,000	£24,000	£48,000	£72,000

TREATMENT	Treatment Course Price	Monthly Courses Treatments	Monthly Revenue Month 1	Month 3	Month 6	Month 9
Full Face Treatment (Single Treatment)	£1,200	1	£1,200	£3,600	£7,200	£10,800
	£1,200	2	£2,400	£7,200	£14,400	£21,600
	£1,200	5	£6,000	£18,000	£36,000	£54,000

TECHNICAL SPECIFICATION

TECHNICAL SPECIFICATION	
Energy Source	RF (Radio Frequency)
Working mode	Bipolar RF
Frequency	1MHz
Output energy	2-120J/cm²
Frequency	1Hz-10Hz
RF Output Power	35W
Pin Number	36 pins (Face) 64 pins (Body)
Pin Diameter	0.12mm
Penetration depth	0.5-3.5mm
Pressure	5 levels - Level 1: 42kpa, Level 2: 46kpa, Level 3: 50kpa, Level 4: 54kpa, Level 5: 60kpa
GENERAL SPECIFICATION	
Power supply	220V-240V, 50/60Hz
Fuse	5 Amp
Max actual power noise	60dB
Treatment time	From 25 to 35 minutes
Visualisation	Graphic Display 10.4" touch screen
Size of machine (h, w, d)	110 x 59 x 54cm
Net weight	42kg
Manufacturer Guarantee	12 months, subject to T&Cs



TRAINING



WHAT DOES THE TRAINING INCLUDE?

- A one day virtual training course.
- A comprehensive training pack with all required training materials.
- The theory behind the technologies.
- Interactive practical sessions led and observed by our 3D trainers.
- VTCT ITEC accredited training*.
- First-hand advice from trainers on using the devices in real clinic settings.
- Unlimited repeat or top-up training sessions as required.
- Guidance on supporting products.



TRAINING



WHAT DO THE TRAINING COURSES COVER?

Each course starts with a morning of theory work. This is where practitioners will learn the science behind the technologies on their machine, in addition to:

- The consultation process.
- Contra-indications and medications.
- 3D approved treatment protocols.
- Different client types and how to treat them.
- Advice on packaging treatments (on our multi-technology devices).

Practical Session:

- Client preparation and consultation.
- Live practical demonstration.
- Student practical application observed by the trainer.
- Pre and post-care information



3D Dermaforce

ADVANCED RADIOFREQUENCY MICRONEEDLING

To find out more or to proceed, get in touch with your dedicated Business Development Manager today.

3D-AESTHETICS.CO.UK

 @3DlipoLtd  @3Dlipo  @3daesthetics_

